

Jazzy® Ultra Light¹

HCPCS Code	K0899 Power Mobility Device, not coded by DME PDAC or does not meet criteria
Weight Capacity	300 lbs.
Turning Radius ³	33.19"
Overall Width ³	24.49"
Overall Length ³	40.91"
Maximum Speed ²	Up to 3.69 mph
Ground Clearance ³	2.17" at motor
Casters Wheels ³	8" Flat free
Drive Wheels ³	10.71" Flat free
Rear Casters ³	N/A
Drivetrain	Two motor, rear-wheel drive
Braking System	"Intelligent Braking" electronic regenerative, disc park brake
Suspension Type	Suspension Caster Forks
Seat Dimensions ³	18" x 18"
Component Weights ⁶	Complete Assembly: 33 lbs. 10Ah Battery: 5 lbs. 20Ah Battery: 7 lbs.
Total Weight ⁴	38 lbs. with 10 Ah Battery 40 lbs. with 20 Ah Battery 33 lbs. without battery
Range Per Charge ^{2,3,5}	10Ah: Up to 8 miles at 200 lbs. Up to 7 miles at 300 lbs. 20Ah: Up to 16.4 miles at 200 lbs. Up to 14.4 miles at 300 lbs.
Battery Requirements	10Ah, 25.2V Lithium-ion 20Ah, 25.2V Lithium-ion
Battery Charger	Off-Board
Warranty	Frame: lifetime limited Drivetrain: 1-year Electronics: 1-year Battery: 13-Month

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.



PRIDE
MOBILITY

(US) 800.800.8586 • (Canada) 888.570.1113

pridemobility.com |      

- ¹ Pride® FDA Class II Medical Devices are designed to aid individuals with mobility impairments.
- ² Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- ³ Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- ⁴ AGM or gel cell type required.
- ⁵ Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- ⁶ Battery weight may vary based on manufacturer.

010825

Jazzy
ULTRA LIGHT



LIGHTWEIGHT



Weights just
33 Lbs
without the
battery

TRAVEL



REDEFINED

pridemobility.com |      

Jazzy®

ULTRA LIGHT

Weighs just
33 Lbs
without the
battery

Breathable
Mesh back

Programmable
joystick controller

Ultra-lightweight,
carbon fiber frame

Standard
mesh cupholder

Airline compliant**
lithium battery
(10Ah Only)

Perforated
Tires

Suspension
caster fork

PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS

FDA Class II Medical Device

JAZZY® KEEPS GETTING LIGHTER!

The **Jazzy® Ultra Light** power wheelchair weighs a remarkable 33 pounds without the battery and offers superior durability thanks to an all-carbon fiber frame. Foldable, portable, and designed for comfort, it's your ideal travel companion for all the journeys that lie ahead.

Features

- Easy fold design that locks into a folded position
- Weighs just 33 pounds without the battery
- Mesh back and seat with standard cushion for extra support
- Durable all-carbon fiber frame
- Programmable joystick controller can attach to either armrest
- Powered by an airline compliant** lithium battery (10Ah Only)
- Ground clearance of 2"
- Front suspension to ensure a smooth ride over various terrain
- Compatible with our new Jazzy Travel Bag
- Rear lights illuminate when unit is in motion
- Mesh cup holder comes standard

Available Color

Black

Accessories

- Jazzy® Travel Bag
- Additional battery
- Headrest and Attendant control bracket bundle
- Padded armrest covers


** Transport of this item on commercial aircraft is subject to federal regulations and airline policy. Please consult your airline prior to flight to ensure your Jazzy® Ultra Light and its lithium-ion battery are permitted.



Perforated Tires



Easy Fold Design that locks into a folded position

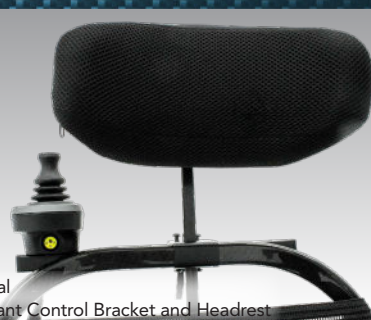

Optional
Attendant Control Bracket and Headrest